

INSPIRED BY NATURE'S  
GOLDEN RATIO

Acriva<sup>LD</sup>  
Trinova  
PRO

pupil adaptive<sup>®</sup>  
Sinusoidal Vision Technology Trifocal IOL

ENHANCED  
SINUSOIDAL TRIFOCAL IOL

# THE GOLDEN RATIO IN ENHANCED SIM



## NO MORE COMPROMISE

Halo & Glare? Not a Problem Anymore<sup>1,2,3</sup>

*Minimized dysphotopsia due to reduced scattered light*



## DISCOVER THE COLORS OF LIFE

Maximum Light Transmission into Retina (93%)<sup>4,5</sup>

*Unique sinusoidal pattern reduces light loss significantly, bringing a higher MTF result and better contrast sensitivity*



# SINUSOIDAL TRIFOCAL IOL TECHNOLOGY



## NO MORE BOUNDARIES

**Comfortable Reading  
Distances<sup>6,7,8,9</sup>**

*Sinusoidal Vision Technology (SVT®)  
provides spectacle-free life from all distances*



## DETAILS ARE IMPORTANT

**Optical Excellence<sup>10</sup>**  
*Outstanding visual outcomes even in  
mesopic conditions*

# ENHANCED SINUSOIDAL TRIFOCAL IOL TECHNOLOGY

Enhanced Sinusoidal Vision Technology developed by VSY Biotechnology is inspired by the **sinusoidal pattern** seen in nature.

## What are Sinusoidal Diffraction Benefits?

The sinusoidal diffraction pattern **has naturally three different foci** appear as  $-1^{\text{st}}$  (works for far),  $0^{\text{th}}$  (works for intermediate), and  $+1^{\text{st}}$  (works for near) diffraction orders (Fig A)<sup>11</sup>. That is why sinusoidal diffraction **does not need to have an overlapping pattern and has a lower ring number**.

- **1 ring for 3 different foci**
- **Less glare and halo, more efficiency<sup>1,2,3</sup>**

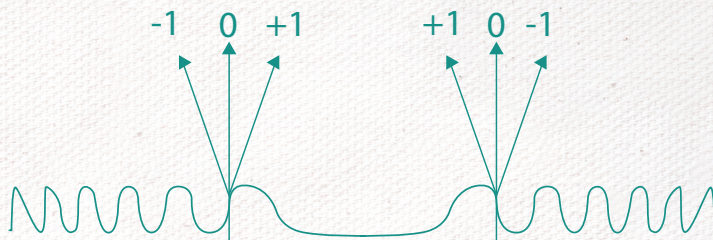


Fig A

No overlapping pattern and 3 different foci  
**Enhanced light utilization<sup>12</sup>**

Enhanced Sinusoidal Vision Technology that mimics a sine wave-like surface profile is designed to obtain the **ideal optical performance**: minimum dysphotopsia, maximum light transmission, and optimum light distribution to let your patients enjoy seamless, continuous vision from all distances day and night.

# NO MORE COMPROMISE

## Halo & Glare? Not a Problem Anymore



The overlapping pattern's sharp edges and pointy peaks on traditional trifocal IOLs are the main causes of positive dysphotopsia.

**Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup> has a smoothly varying surface profile that helps to reduce halos and glare by reducing scattered light.**

**Acriva<sup>UD</sup>**  
**Trinova**  
**PRO**  
pupil adaptive<sup>®</sup>  
Sinusoidal Vision Technology Trifocal IOL



**Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>,  
Sinusoidal Trifocal IOL Surface Profile  
Unique Sinusoidal Pattern**

## TRADITIONAL TRIFOCAL IOL




**Traditional Trifocal IOL Surface Profile  
Overlapping Pattern with Sharp Edges**

# DISCOVER THE COLORS OF LIFE

## Maximum Light Transmission into Retina (93%)

*Our unique sinusoidal pattern reduces light loss significantly, bringing a higher MTF result and better contrast sensitivity.*

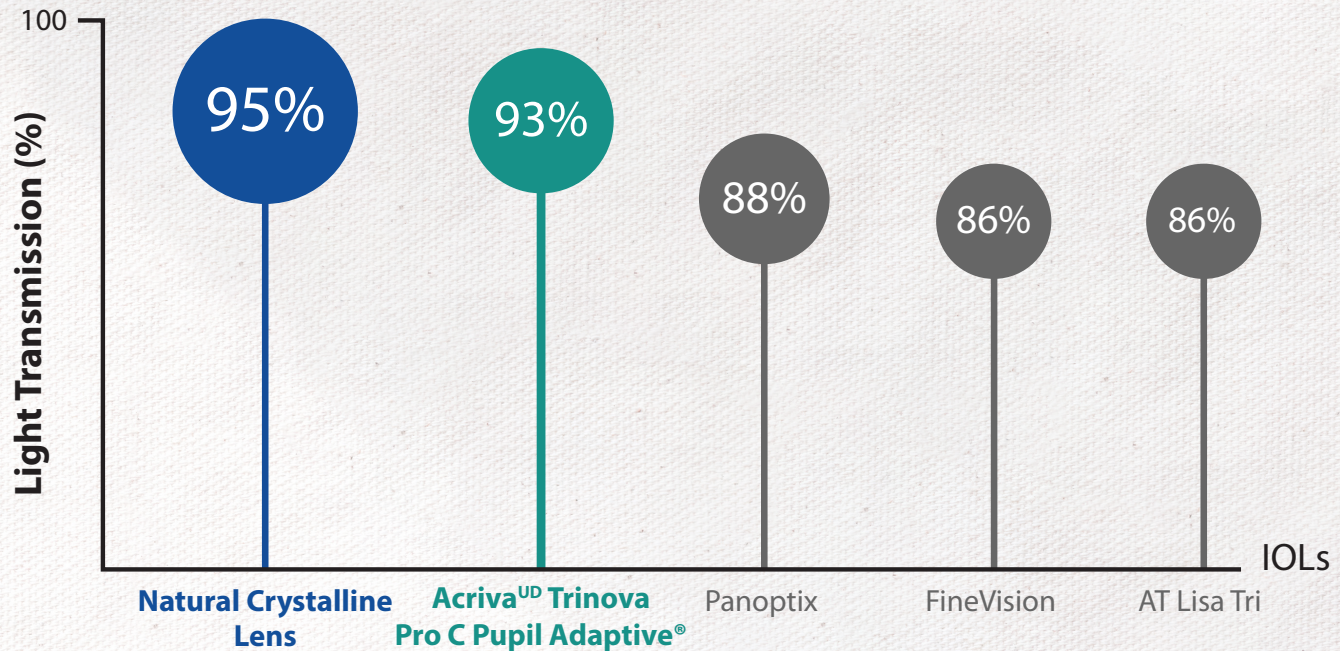


The Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>'s Enhanced Sinusoidal Optical profile is designed to ensure maximum light transmission (93%) to the retina, providing improved contrast sensitivity to enable your patients to enjoy all the colors of life.

The crystalline lens of a healthy 30-year-old individual has a 95% light transmission and Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup> offers the closest light transmission to this rate (93%). Trifocal IOLs with blazed diffractive patterns are known to cause significant light loss. Even 1% of light loss affects a patient's overall visual performance exponentially.

# NO MORE BOUNDARIES

## Comfortable Reading Distances

Enhanced Sinusoidal Vision Technology provides spectacle-free life from all distances.

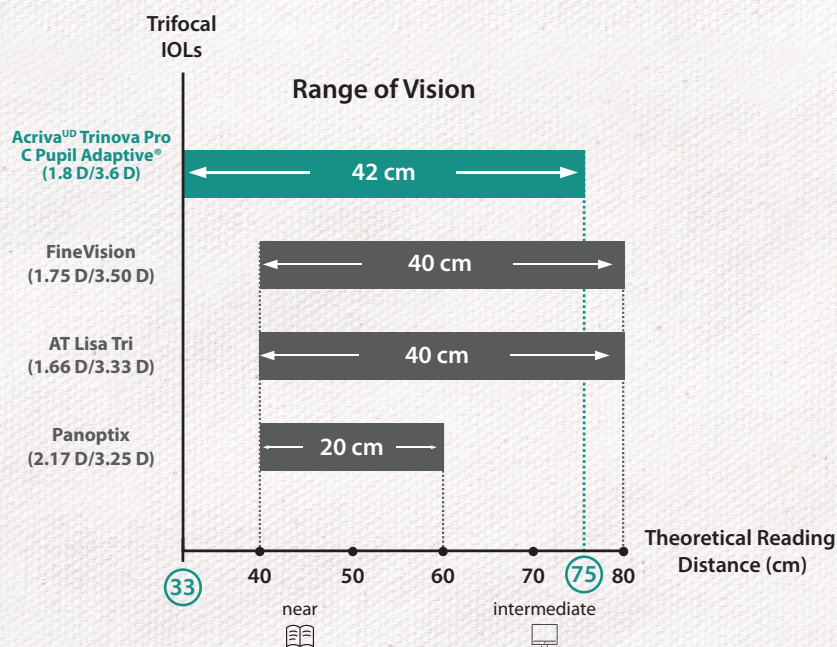

**+3.6 D** near and **+1.8 D** intermediate addition of Acryva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup> design offers exceptional contrast sensitivity and visual acuity from near, intermediate, and far distances.

Enhanced Sinusoidal Vision Technology is designed to provide the seamless vision required for all daily activities (reading, using mobile phone, watching TV, cooking, driving, etc.) without spectacles.



Acryva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup> offers the widest range of vision (between 33- 75 cm) compared to other trifocal IOLs listed in the table.

# DETAILS ARE

## Optical Exo

Outstanding visual outcomes e

**Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>** is designed to produce the best possible MTF results in a trifocal IOL. The MTF outcome shows the image quality at different focal points.

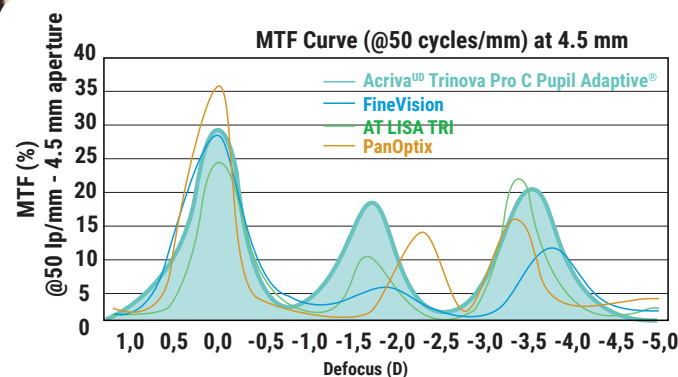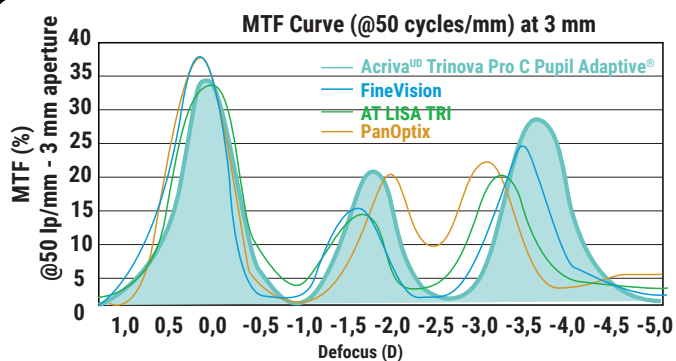

As seen in the MTF curve, at photopic conditions, **Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>** showed exceptional performance throughout the foci thanks to its unique sinusoidal optic design.

In light of the MTF curve in mesopic conditions, **Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>** clearly provides better optical image quality—namely, better contrast at near and intermediate distances compared to other trifocal IOLs.



# IMPORTANT

excellence

even in mesopic conditions

|  | 3 mm |              |      | 4.5 mm |              |      |
|--|------|--------------|------|--------|--------------|------|
|  | FAR  | INTERMEDIATE | NEAR | FAR    | INTERMEDIATE | NEAR |
| <b>Acriva<sup>®</sup><br/>Trinova<br/>PRO</b><br>pupil adaptive <sup>™</sup><br><small>Simultaneous Vision Technology • Trifocal IOL</small> |      |              |      |        |              |      |
| FineVision   |      |              |      |        |              |      |
| AT LISA Tri  |      |              |      |        |              |      |
| PanOptix   |      |              |      |        |              |      |

As the name implies, a resolution test pattern is a tool for measuring the resolving power of an optical system. It consists of reference line patterns with well-defined thicknesses and spacings. The test helps demonstrate the performance of a lens's resolution in photopic and mesopic conditions. The USAF resolution test performed for **Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>** at different apertures revealed excellent performance compared with other trifocal IOLs.



# GOLDEN RATIO

## ENERGY DISTRIBUTION A

Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup> smartly distributes the light energy across different pupil sizes to maximize light utility, retinal photoreceptor interactions, high contrast, and visual distance, eliminating the need for spectacles.

### NEAR

### INTERMEDIATE

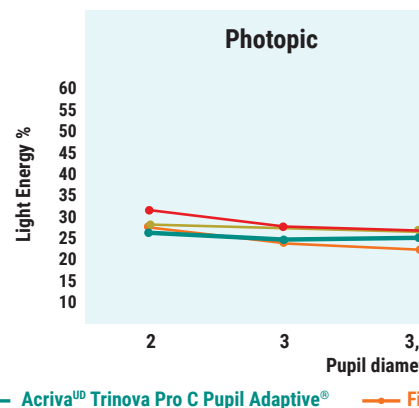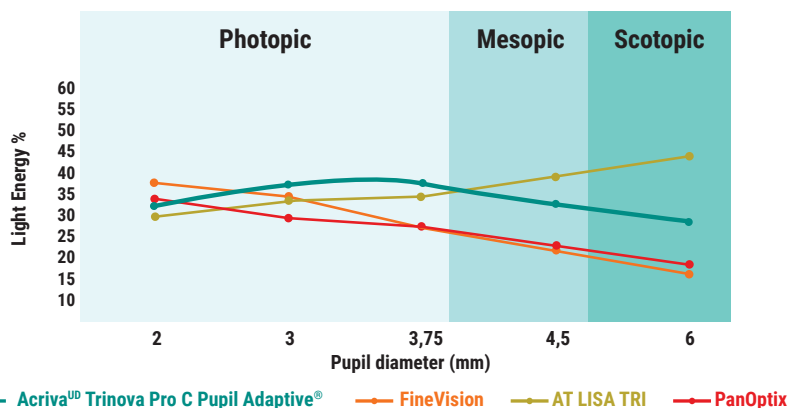

Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup> distributes the light energy between the foci in such a way so that in mesopic illumination, a significant amount of light energy is distributed to far and intermediate vision more than near, as shown in the above figures in mesopic conditions.

Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup> distributes effective light energy across different foci while maintaining distance visual acuity due to minimal light loss and enhanced light distribution with patented technology.




# IN TRIFOCALITY

## ADAPTED BY PUPIL SIZE


energy between the three foci according to different illuminance levels and high contrast sensitivity, resolution and wide range of vision including reading spectacles in every lighting situation.

### MEDIATE

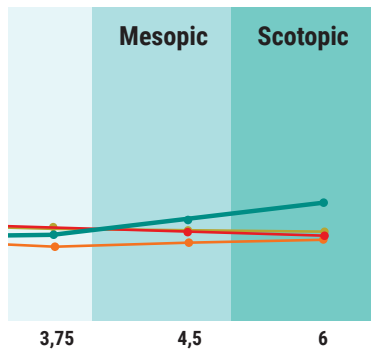

3,75 4,5 6  
diameter (mm)  
— FineVision — AT LISA TRI — PanOptix

### FAR

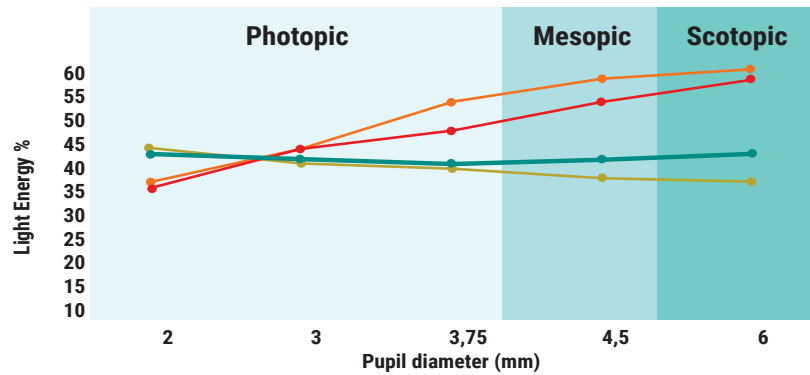

— Acryva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup> — FineVision — AT LISA TRI — PanOptix

illuminance levels, the light energy is directed to the near and intermediate vision, while in scotopic illuminance levels, this mesopic and scotopic conditions, and near vision is empowered for photopic conditions.

ive and efficient light energy into the near and intermediate  
minimized light loss in the diffraction orders provided by  
tentented Sinusoidal Vision Technology.



# GOLDEN RATIO IN TRIFOCALITY: IOL ADAPTING TO REAL LIFE CONDITIONS...


**Pupil diameter changes by lighting condition and eye accommodation state.**<sup>13,14,15</sup>

**Changing pupil diameter affects our visual performance due to the retinal illumination level and the directional sensitivity of the retinal photoreceptor cells.**<sup>16</sup>

Retinal image quality has been shown to correlate better with VA tasks of lower contrast and lower luminance than the task of photopic high contrast (HC)-VA.

Therefore evaluating IOL performance and patient satisfaction in a standard examination room under the ideal conditions of 2-3 mm pupil diameter, mega contrast letters, and optimum light conditions do not reflect the real-life experience.<sup>17,18,19,20</sup>

**Now Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>, with enhanced light distribution with patented Sinusoidal Vision Technology, each sinusoidal wave is optimized to reach an ideal design just like the Golden Ratio in nature. Thanks to renewed unique design, the lens adapts to the pupil diameter and provides maximum light efficiency, optimum energy distribution, and maximum visual acuity.**

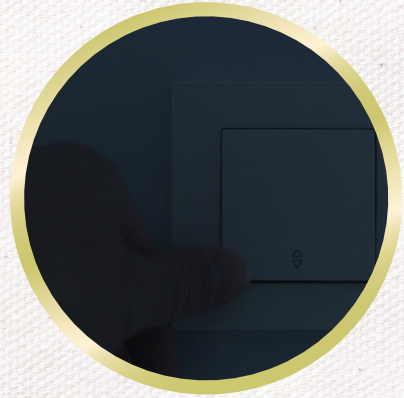

**Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>** is designed for today's lifestyles; from viewing traffic signs comfortably and clearly while driving at night safely, to using mobile devices and computer screens to reading a book comfortably even in dim light with high-quality distance vision in a range of different light conditions. Thus it maximized spectacle independence at all distances and in every possible lighting condition.

**Acriva<sup>UD</sup> Trinova Pro C Pupil Adaptive<sup>®</sup>** aims to improve patients' quality of life by maximizing the quality of vision.

# TECHNICAL


 **Technology**  
Enhanced Sinusoidal  
Pupil Adaptive Trifocal IOL

 **Optic Size**  
6.00 mm

 **Overall Length**  
13.00 mm

 **Haptic Design**  
C Loop (suitable for MICS)

 **Haptic Angle**  
0°


 **Material**  
Hydrophilic acrylic with 25% water  
content, Natural Chromophore with UV,  
Violet, Blue Light Filter

 **Asphericity**  
Mildly Negative Bi-Aspheric Design

 **Abbe Number**  
58

# FEATURES

 **Light Transmission**  
93%

 **Dioptric Power Range**  
Sph. 0.0 to +32.0 D (0.5 D increments)

 **Square Edge Design**  
360° All Enhanced Square Edge

 **Nominal A-Constant**  
118.0

 **Refractive Index**  
1.46

## References

1. Moreno, V., Román, J. F., & Salgueiro, J. R. (1997). High efficiency diffractive lenses: Deduction of kinoform profile. *American Journal of Physics*, 65(6), 556–562.
2. Stodulka P. Clinical results of implantation of the new sinusoidal trifocal iol . Presented at the 36th congress of ESCRS symposium, Vienna Austria September 2018
3. Tomita, M. (2014). Diffractive Multifocal IOLs: The Acryva<sup>®</sup> Reviol MFM 611 IOL and Acryva<sup>®</sup> Reviol MF 613 IOL. In *Multifocal Intraocular Lenses* (pp. 147-153). Springer, Cham.
4. Vega, F., Valentino, M., Rigato, F., & Millán, M. S. (2021). Optical design and performance of a trifocal sinusoidal diffractive intraocular lens. *Biomedical Optics Express*, 12(6), 3338-3351.
5. Albero, J., Davis, J. A., Cottrell, D. M., Granger, C. E., McCormick, K. R., & Moreno, I. (2013). Generalized diffractive optical elements with asymmetric harmonic response and phase control. *Applied optics*, 52(15), 3637-3644.
6. Ceran, B. B., Arifoglu, H. B., Ozates, S., & Tasindi, E. E. (2020). Refractive results, visual quality and patient satisfaction with a new trifocal intraocular lens design. *Annals of Medical Research*, 27(11), 3018-3023.
7. Mrukwa-Kominek, E., et al. "The sinusoidal trifocal intraocular lens in cataract surgery and its effect on the quality of patients' vision." ESCRS, Paris 2019
8. Kontadakis, G., et al. "Visual acuity and patient satisfaction after bilateral implantation of a trifocal enhanced-depth-of-focus intraocular lens". ESCRS, Paris, 2019
9. Europe, Multicentric trials, VSY Biotechnology Data on File (2021)
10. VSY Biotechnology R&D Center, Data on File (2021)
11. Sokolowski, M., Pniewski, J., Brygola, R., & Kowalczyk-Hernandez, M. (2015). Hybrid hepta-focal intraocular lenses. *Optica Applicata*, 45(3).
12. Valle, P. J., Oti, J. E., Canales, V. F., & Cagigal, M. P. (2005). Visual axial PSF of diffractive trifocal lenses. *Optics express*, 13(7), 2782-2792.
13. Fry, G. A. (1945). The relation of pupil size to accommodation and convergence. *Optometry and Vision Science*, 22(10), 451-465.
14. Napieralski, P., & Rynkiewicz, F. (2019). Modeling human pupil dilation to decouple the pupillary light reflex. *Open Physics*, 17(1), 458-467.
15. Plakitsi, A., & Charman, W. N. (1997). Ocular spherical aberration and theoretical through-focus modulation transfer functions calculated for eyes fitted with two types of varifocal presbyopic contact lens. *Contact Lens and Anterior Eye*, 20(3), 97-106.
16. Westheimer, G. (2008). Directional sensitivity of the retina: 75 years of Stiles–Crawford effect. *Proceedings of the Royal Society B: Biological Sciences*, 275(1653), 2777-2786.
17. Puell, M. C., Perez-Carrasco, M. J., Palomo-Alvarez, C., Antona, B., & Barrio, A. (2014). Relationship between halo size and forward light scatter. *British Journal of Ophthalmology*, 98(10), 1389-1392.
18. Ravalico, G., Baccara, F., & Rinaldi, G. (1993). Contrast sensitivity in multifocal intraocular lenses. *Journal of Cataract & Refractive Surgery*, 19(1), 22-25.
19. Tanabe, H., Tabuchi, H., Shoji, T., Yamauchi, T., & Takase, K. (2020). Comparison of visual performance between monofocal and multifocal intraocular lenses of the same material and basic design. *Scientific reports*, 10(1), 1-11.
20. Das, K. K., Stover, J. C., Schwiegerling, J., & Karakelle, M. (2013). Technique for measuring forward light scatter in intraocular lenses. *Journal of Cataract & Refractive Surgery*, 39(5), 770-778.

**Acriva**<sup>LD</sup>  
**Trinova**  
**PRO**  
pupil adaptive<sup>®</sup>  
Sinusoidal Vision Technology Trifocal IOL



**VSY Biotechnology GmbH**  
Esslinger Str. 7, 70771 Leinfelden-Echterdingen, Germany  
contact@vsybiotechnology.com / www.vsybiotechnology.com

